



[My Account](#) | [Order Status](#) | [Quick Order](#) [Worldwide](#)



- [Products & Services](#)
 - [Applications](#)
 - [Nucleic Acid Purification & Analysis](#)
 - [Nucleic Acid Amplification & Expression Profiling](#)
 - [Cloning](#)
 - [Proteins, Expression, Isolation & Analysis](#)
 - [RNAi, Epigenetics & Gene Regulation](#)
 - [Cell Culture](#)
 - [Cell & Tissue Analysis](#)
 - [Drug Discovery](#)
 - [Clinical & Diagnostic Applications](#)
 - [Bioproduction](#)
 -
 - [Product Types](#)
 - [Antibodies & Secondary Detection Reagents](#)
 - [Beads & Microspheres](#)
 - [Buffers & Chemicals](#)
 - [Cell Culture & Transfection Reagents](#)
 - [Columns, Resins & Plates](#)
 - [Competent Cells & Strains](#)
 - [Controls, Standards & Ladders](#)
 - [DNA Vectors, Clones & Libraries](#)
 - [Eukaryotic Cells](#)
 - [Gels, Fractionation Strips, & Membranes](#)
 - [Instruments, Software, & Accessories](#)
 - [Kits & Assays](#)
 - [Labeling & Detection Products](#)
 - [Literature](#)
 - [Primers, Oligos, & Nucleotides](#)
 - [Proteins, Enzymes, & Peptides](#)
 - [Slide & Chip Arrays](#)
 -
 - [Services](#)
 - [Discovery Research](#)
 - [Production & Process Development](#)
 - [OEM Services](#)
 - [Cell Culture Services™](#)
 - [Licenses, Terms, and Conditions](#)
 -
 - [Promotions](#)
 - [New Products](#)
 - [Sample Offers](#)
 - [Check Price Tool](#)
 - [Quick Order Tool](#)
 - [Supply Center](#)
 - [Request Information](#)

- [SCMS Account Registration](#)
 - [SCMS Admin Login](#)
 -
 - [eProcurement](#)
 -
-
- [Applications](#)
 - [Nucleic Acid Purification & Analysis](#)
 - [DNA Purification](#)
 - [RNA Purification](#)
 - [Viral RNA/DNA Purification](#)
 - [Nucleic Acid Quantitation](#)
 - [Nucleic Acid Gel Electrophoresis](#)
 - [Automated Nucleic Acid Purification](#)
 -
 - [Nucleic Acid Amplification & Expression Profiling](#)
 - [Oligonucleotide Design](#)
 - [Real-Time PCR](#)
 - [PCR](#)
 - [Genomic Profiling - Array Labeling & Genotyping](#)
 - [Reverse Transcription & cDNA Synthesis](#)
 - [Gene Expression Profiling - Array Labeling & SAGE](#)
 - [RT-PCR](#)
 - [DNA Sequencing](#)
 - [Nucleic Acid Labeling](#)
 -
 - [Cloning](#)
 - [PCR Cloning](#)
 - [PCR](#)
 - [Gateway® Cloning](#)
 - [Transformation](#)
 - [Restriction & Cloning Enzymes](#)
 - [Oligonucleotide Design](#)
 - [cDNA Libraries & Library Construction](#)
 - [Cloning Services](#)
 - [Mutagenesis](#)
 - [Premade Clones](#)
 - [Microbiological Culture](#)
 - [Vector NTI® Software](#)
 -
 - [Proteins, Expression, Isolation & Analysis](#)
 - [Antibodies and Purified Proteins](#)
 - [ELISA](#)
 - [Luminex® Assays](#)
 - [Protein Expression Systems & Vectors](#)
 - [Protein Quantitation](#)
 - [Protein-Protein Interactions](#)
 - [Protein Sample Preparation & Protein Isolation](#)
 - [Mass Spectrometry](#)
 - [Enzyme & Protein Activity Assays](#)
 - [Transfection & Selection](#)
 - [Protein Gel Electrophoresis](#)
 - [Western Blotting](#)


- [Biomarker Discovery](#)
 -
- [RNAi, Epigenetics & Gene Regulation](#)
 - [RNAi](#)
 - [miRNA Profiling](#)
 - [Methylation Analysis](#)
 - [Transcription Factor Binding Assays](#)
 - [Chromatin Immunoprecipitation \(ChIP\)](#)
 -
- [Cell Culture](#)
 - [Mammalian Cell Culture](#)
 - [Stem Cell Research](#)
 - [Transfection & Selection](#)
 - [Insect Cell Culture](#)
 - [3D Cell Culture](#)
 - [Cell Counting, Viability & Cryopreservation](#)
 - [Cell Isolation & Expansion](#)
 - [Primary Cell Culture](#)
 -
- [Cell & Tissue Analysis](#)
 - [Antibodies & Secondary Detection](#)
 - [Cell Tracing, Tracking](#)
 - [Cell Isolation & Expansion](#)
 - [Cellular Imaging](#)
 - [Stem Cell Differentiation](#)
 - [Flow Cytometry](#)
 - [Cell Viability, Proliferation & Function](#)
 - [Immunoassays](#)
 - [Cell Structure](#)
 - [Labeling Chemistry](#)
 - [Quantum Dots & Microspheres](#)
 - [Fluorescent Microspheres](#)
 - [Microbiological Analysis](#)
 -
- [Drug Discovery](#)
 - [Target & Lead Identification & Validation](#)
 - [Drug Metabolism & Safety](#)
 -
- [Clinical & Diagnostic Applications](#)
 - [Transplant Diagnostics](#)
 - [Clinical Cell Isolation & Expansion](#)
 - [Anatomical Pathology](#)
 - [Cytogenetics](#)
 - [BioDefense & Pathogen Detection](#)
 - [Environmental Diagnostics](#)
 - [Bead-Based IVD Assays](#)
 -
- [BioProduction](#)
 - [Protein & Antibody Production](#)
 - [Virus & Vaccine Production](#)
 - [Custom Media Products](#)
 - [Bioproduction Services](#)
 - [Basic Cell Culture](#)

-
-
-
- Brands
 - Dynal®
 - Dynabeads® Products & Technology
 - The History of Dynal® and Biomagnetic Separation
 - Bead Types
 - Magnets
 -
 - Gibco
 - New Technologies
 - General Information
 - Reference Guide
 - Press Releases
 -
 - Molecular Probes®
 - The Handbook
 - Literature
 - Subscribe to ProbesOnline Newsletter
 - Legal Notices and Trademark Attributions
 - Licensing and Business Development
 -
 - Product Brands
 -
-
- Guides
 - LINNEA™ Communities
 - Vector NTI® Community
 -
 - LINNEA™ Designers
 - BLOCK-iT™ RNAi Designer
 - D-Lux™ Designer
 - OligoPerfect™ Designer
 - PeptideSelect™
 - VectorDesigner™
 -
 - LINNEA™ Genes
 - LINNEA™ Guide to Antibodies
 - LINNEA™ Guide to Clones
 - Ultimate™ ORF Clones
 - Clone Search
 -
 - LINNEA™ Guide to Primers
 - D-Lux™ Assay Search
 - D-Lux™ Designer
 - OligoPerfect™ Designer
 -
 - LINNEA™ Guide to RNAi
 - BLOCK-iT™ RNAi Express
 - BLOCK-iT™ RNAi Designer
 -
 - LINNEA™ Pathways

- [LINNEA™ Protocols](#)
-
-
- [Support](#)
 - [Order Support - FAQ](#)
 - [Product Support - FAQ](#)
 - [Scientific Tools & Tips](#)
 - [Newsletters & Journals](#)
 - [Cell Culture Archive](#)
 - [Benchtopics Archive](#)
 - [Focus Journal Archive](#)
 - [Expressions Archive](#)
 - [Postings Archive](#)
 - [Quest Archive](#)
 -
 - [Research Tools](#)
 - [KD calculator](#)
 - [Image Gallery](#)
 - [Base-Dye Ratio Calculator](#)
 - [Fluorescence SpectraViewer](#)
 -
 - [Online Seminars](#)
 - [Product Technical Resources](#)
 - [Product Spectra](#)
 - [Product Structures](#)
 - [Citations and Bibliographic References Search](#)
 - [Cell Lines Advanced Search](#)
 - [Cell Lines Browse](#)
 -
 - [Tutorials](#)
 - [Search Tips](#)
 -
-
- [Corporate](#)
 - [Who We Are](#)
 - [Our Quest](#)
 - [Company Fact Sheet](#)
 - [Areas of Focus](#)
 - [Global Locations & Directions](#)
 - [Awards & Recognition](#)
 - [Customer Testimonials](#)
 -
 - [Corporate Citizenship](#)
 - [Community Involvement](#)
 - [Invitrogen Foundation](#)
 - [Environmental Responsibility](#)
 -
 - [Careers](#)
 - [Job Listings](#)
 - [Inside Invitrogen](#)
 - [Meet Our Employees](#)
 - [Benefits](#)
 -

- [Media Center](#)
 - [Press Releases](#)
 - [Image Gallery](#)
 - [Videos & Podcasts](#)
 - [Media Contacts](#)
 - [News Subscription](#)
 -
- [Investor Relations](#)
 - [All SEC Filings](#)
 - [Calendar of Events](#)
 - [Corporate Governance](#)
 - [Investor FAQs](#)
 - [Financials](#)
 - [Investor Kits](#)
 - [Presentations](#)
 - [Press Releases](#)
 - [Shareholder Information](#)
 - [Stock Information](#)
 - [Web Cast Archive](#)
 -
-

Search

All Results 

[Search](#)

View Cart  0 Items

LOGIN 

Molecular Probes

- [The Handbook](#)
 - [Recent Updates](#)
 - [Technical Notes and Product Highlights](#)
 - [List of Tables](#)
 - [List of Figures](#)
- [Literature](#)
- [Subscribe to ProbesOnline Newsletter](#)
- [Subscribe to Our Mailing List](#)
- [Legal Notices and Trademark Attributions](#)
- [Licensing and Business Development](#)

Have a question or
just need help?

Contact us for
Service and Support




[Home](#) > [Brands](#) > [Molecular Probes](#) > [Literature](#) > The Handbook

Chapter Section Section Chapter


Chapter 8 — Nucleic Acid Detection and Genomics Technology

Section 8.1 — Nucleic Acid Stains

 [Printer-Friendly Format](#)

- Properties of Cyanine Dyes
- Premier Cyanine Dyes for Ultrasensitive Nucleic Acid Detection and Quantitation
- Cell-Impermeant Cyanine Dimers: The TOTO Family of Dyes
 - High Affinity for Nucleic Acids
 - High Fluorescence Enhancements and High Quantum Yields upon Binding to Nucleic Acids
 - Modifying the Dimers Creates Compounds with Different Spectral Characteristics
 - Binding Modes of the Cyanine Dimers
 - Working with Cyanine Dimers
- Cell-Impermeant Cyanine Monomers: The TO-PRO Family of Dyes
 - Spectral Characteristics of the Cyanine Dye Monomers
 - Working with Cyanine Monomers
- Cell-Impermeant SYTOX Dyes for Dead-Cell Staining
 - SYTOX Green Stain
 - SYTOX Blue Stain
 - SYTOX Orange Stain
- Cell-Permeant Cyanine Dyes: The SYTO Nucleic Acid Stains
 - SYTO Nucleic Acid Stains for DNA and RNA
 - SYTO RNASelect Green-Fluorescent Cell Stain
- Chemically Reactive Cyanine Dyes
- RediPlate 96 Nucleic Acid Stain Sampler Microplate
- Phenanthridines and Acridines: Classic Intercalating Dyes
 - Cell-Impermeant Ethidium Bromide and Propidium Iodide
 - Cell-Permeant Hexidium Iodide
 - Cell-Permeant Dihydroethidium (Hydroethidine)
 - High-Affinity Ethidium Homodimers
 - Ethidium Monoazide: A Photocrosslinking Reagent
 - Acridine Orange: A Dual-Fluorescence Nucleic Acid Stain
 - AT-Selective Acridine Homodimer
 - AT-Selective ACMA
- Indoles and Imidazoles: Classic Minor Groove-Binding Dyes
 - DNA-Selective Hoechst Dyes
 - AT-Selective DAPI
- Other Nucleic Acid Stains
 - 7-Aminoactinomycin D and Actinomycin D: Fluorescent Intercalators
 - Multicolor Hydroxystilbamidine
 - Long-Wavelength LDS 751
 - NeuroTrace Fluorescent Nissl Stains
- [Data Table](#)
- [Product List](#)


Section 8.2 — Labeling Oligonucleotides and Nucleic Acids

 [Printer-Friendly Format](#)

- ChromaTide Nucleotides
 - Structures of the ChromaTide Nucleotides
 - Fluorescent ChromaTide Nucleotides
 - ChromaTide Dinitrophenyl (DNP)-11-dUTP
 - Using ChromaTide Nucleotides in Enzymatic Labeling Methods
- Amine-Modified Nucleotides
 - Unlabeled and Labeled aha-dUTP and aha-dCTP
 - Aminoallyl UTP and Aminoallyl dUTP
 - Alexa Fluor Amine-Reactive Dye Decapacks for Labeling Amine-Modified DNA and RNA


- 5-Bromo-2'-Deoxyuridine, 5-Bromo-dUTP (BrdUTP) and 5-Bromo-UTP (BrUTP)
- ULYSIS Nucleic Acid Labeling Kits
- ARES DNA Labeling Kits
- Labeled Oligonucleotides
- Labeling Amine- and Thiol-Modified Oligonucleotides
 - Alexa Fluor Oligonucleotide Amine Labeling Kits
 - Alexa Fluor Amine-Reactive Dye Decapacks
 - Other Reactive Fluorescent Dyes and Quenchers
 - Dyes for Sequencing Applications
- Labeling Phosphate-Modified Oligonucleotides
- Other Chemical Labeling Methods for Nucleic Acids
 - Labeling Cytidine Residues
 - Labeling the 3'-Terminus of RNA
 - Labeling Abasic Sites with ARP
 - Specialized Methods for Nucleic Acid Modification
- ElutaTube Microdialysis Vials
- [Product List](#)

Section 8.3 — Nucleic Acid Detection and Quantitation in Solution

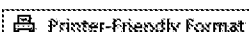


- Nucleic Acid Detection and Quantitation in Solution
- Quant-iT Assay Kits for DNA and RNA
 - Quant-iT DNA Assay Kits
 - Quant-iT RNA Assay Kit
- PicoGreen dsDNA Quantitation Assay
 - PicoGreen dsDNA Quantitation Reagent and Kits
 - RediPlate 96 PicoGreen dsDNA Quantitation Kits and Microplate
- OliGreen ssDNA Quantitation Reagent and Kit
- RiboGreen RNA Quantitation Assay
 - RiboGreen RNA Quantitation Reagent and Kit
 - RediPlate 96 RiboGreen RNA Quantitation Kit
 - RiboGreen RNA-Specific Quantitation Kit with DNase I
- Other Stains for Nucleic Acid Quantitation in Solution
 - Cyanine Dyes and Phenanthridine Dyes for Nucleic Acid Quantitation in Solution
 - Hoechst 33258 Dye for Quantitating DNA in Solution
- Real-Time Quantitative PCR Using the SYBR Green I Nucleic Acid Gel Stain
- A Useful Technique to Reduce Background Fluorescence in Nucleic Acid Detection
- [Product List](#)

Section 8.4 — Nucleic Acid Detection and Quantitation in Electrophoretic Gels and Capillaries

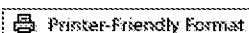

- Nucleic Acid Detection in Gels
 - SYBR Gold Nucleic Acid Gel Stain
 - SYBR Green I Nucleic Acid Gel Stain
 - SYBR Green II RNA Gel Stain
 - SYBR Green Nucleic Acid Gel Stains: Special Packaging and a Starter Kit
 - SYBR Safe DNA Gel Stain
 - Ethidium Bromide
 - Cyanine Monomers for Staining DNA in Electrophoretic Gels
 - Cyanine and Ethidium Dimers for Staining DNA Prior to Electrophoresis
 - Electrophoretic Mobility-Shift (Bandshift) Assay (EMSA) Kit
 - Other Nucleic Acid Stains for Gel-Staining Applications
- Safe Imager Blue-Light Transilluminator
- Capillary Electrophoresis and Channel Electrophoresis
 - Capillary Electrophoresis
- [Product List](#)


Section 8.5 — Detecting Nucleic Acid Hybridization

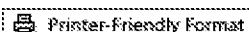

- Principles of Fluorescence *In Situ* Hybridization (FISH)
- Probe Preparation for FISH


- Enzymatic Incorporation of ChromaTide Labeled Nucleotides
 - Aminoallyl dUTP Labeling Using the ARES Kits
 - Direct Chemical Labeling of Nucleic Acids Using the ULYSIS Kits
 - Labeling Oligonucleotides
- Signal Amplification for FISH Using Secondary Detection Reagents
 - Biotin and Streptavidin
 - Signal Amplification Using Anti-Dye Antibodies
 - Tyramide Signal Amplification (TSA)
 - Signal Amplification Using Enzyme-Labeled Fluorescence (ELF) Technology
 - Signal Amplification Using Two Rounds of Amplification
 - Colorimetric Signal Amplification
- Counterstaining Chromosomes for FISH
- Gaskets for *In Situ* Hybridization Experiments
- Nucleic Acid Hybridization on Blots and Microarrays
- Labeling Nucleic Acid Samples for Microarray Experiments
 - Consistent, Uniform Sample Labeling with ARES DNA Labeling Kits
 - Alexa Fluor Reactive Dye Decapacks for Microarray Applications
 - Fast and Easy Direct Labeling with ULYSIS Nucleic Acid Labeling Kits
 - Fluorophore-Labeled Nucleotides
 - Labeled dT₁₈ and Random Oligodeoxynucleotides
 - Using Labeled Oligonucleotides as Biosensors on Solid Supports
- Secondary Detection for Signal Amplification
 - Detecting Biotin Labels with Fluorescent Streptavidin Conjugates
 - Signal Amplification with the ELF 97 Phosphatase Substrate
 - Chemiluminescent Detection on Blots
 - Other Types of Signal Amplification
- Tools and Techniques for Quality Control Testing and Normalization on Microarrays and Blots
 - PARAGON DNA Microarray QC Stain Kit
 - PARAGON DNA Microarray QC Hybridization Kits
 - PARAGON Genomic DNA Hybridization Test Kits
 - PARAGON Microarray Scanner Calibration Slide
 - PARAGON Dye Ratio Calibration Slide
 - Other Nucleic Acid Stains for Standardizing Microarrays
 - SYBR DX DNA Blot Stain
 - Panomer Random-Sequence Oligonucleotides
- Sample Quantitation
- Hybridization Chambers for Microarrays
- Solution-Based Hybridization Assays
 - Real-Time Quantitative PCR Using FRET and Quenching Techniques
 - Other Solution-Based Hybridization Assays Using FRET and Quenching Techniques
 - Other Solution-Based Hybridization Assays Using Fluorescent Nucleic Acid Stains
- [Product List](#)

Section 8.6 — Nuclear and Chromosome Counterstaining and Nissl Stains

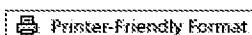

- Nuclear Counterstaining of Fixed Cells and Tissues
 - Blue-Fluorescent Counterstains
 - Green-Fluorescent Counterstains
 - Yellow-Fluorescent Counterstain
 - Orange- and Red-Fluorescent Counterstains
 - Long-Wavelength Nuclear Counterstains
 - SelectFX Nuclear Labeling Kit
 - Cytological Nuclear Counterstain Kit
- Nuclear Counterstains for Live Cells and Unfixed Tissues
 - Cell-Permeant Blue-Fluorescent Counterstains
 - Cell-Permeant Green-Fluorescent Counterstains
 - Cell-Permeant Orange- and Red-Fluorescent Counterstains
- Tracking Chromosomes through Mitosis
- Chromosome Counterstaining
 - Red-Fluorescent Counterstaining for Green-Fluorescent Probes
 - The Best Green-Fluorescent Counterstains
 - Blue-Fluorescent Counterstaining for Multicolor Labeling
 - Very Long-Wavelength Counterstains for Multicolor Labeling
- Chromosome Banding
 - SYTOX Green Nucleic Acid Stain
 - Acridine Homodimer
 - Other Dyes and Chromosome Banding Reagents
- NeuroTrace Fluorescent Nissl Stains
- [Product List](#)

Section 8.7 — Analysis of DNA Structure, DNA Binding and DNA Damage

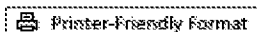

- Nucleic Acid Conformational Analysis
- Examining the Behavior of Single Nucleic Acid Molecules
- DNA Binding Assays
 - Electrophoretic Mobility-Shift (Bandshift) Assays
 - DNA Binding Assays in Solution
 - Selective Cleavage of Nucleic Acids with a Chemical Nuclease
- Assessing DNA Damage
 - Comet (Single-Cell Gel Electrophoresis) Assay to Detect Damaged DNA
 - TUNEL Assay for *In Situ* Detection of Fragmented DNA
 - Microplate Assays for DNA Damage
- Assays for Enzymes that Modify Nucleic Acids
 - Gel Assays for DNase Detection
 - Solution Assays for Nuclease Detection
 - An Assay for Reverse Transcriptase Activity
 - Telomerase
- [Product List](#)

List of Tables

[Table 8.1 — Specialty nucleic acid reagents for molecular biology](#)

[Table 8.2 — Cell membrane-impermeant cyanine nucleic acid stains](#)

[Table 8.3 — Cell-permeant cyanine nucleic acid stains](#)

[Table 8.4 — Properties of classic nucleic acid stains](#)

[Table 8.5 — General characteristics of the dyes provided in the RediPlate 96 nucleic acid stain sampler microplate \(R32715\)](#)

[Table 8.6 — Characteristics of ChromaTide UTP nucleotides](#)

[Table 8.7 — Characteristics of ChromaTide dUTP, ChromaTide OBEA-dCTP, aha-dUTP and aha-dCTP labeled nucleotides](#)

[Table 8.8 — Spectral characteristics of the fluorescent dyes available in Molecular Probes' ULYSIS Nucleic Acid Labeling Kits](#)

[Table 8.9 — Spectral characteristics of the fluorescent dyes available in Molecular Probes' ARES DNA Labeling Kits](#)

[Table 8.10 — Oligonucleotide Amine Labeling Kits](#)

[Table 8.11 — Amine-reactive dyes for nucleic acid sequencing](#)

[Table 8.12 — Selection Guide for the Quant-iT Assay Kits](#)

[Table 8.13 — Effects of Contaminants on the PicoGreen Assay](#)

[Table 8.14 — Effects of Contaminants on the OliGreen Assay](#)

[Table 8.15 — Comparison of the RiboGreen RNA quantitation reagent and ethidium bromide](#)

[Table 8.16 — Effects of Contaminants on the RiboGreen Assay](#)

[Table 8.17 — Tools for hybridization experiments](#)

[Table 8.18 — Spectral characteristics of labeled oligonucleotides](#)

[Table 8.19 — Nucleic acid stains used for microarray quality control](#)

List of Figures

[Figure 8.1 — Fluorescence emission spectra of DNA-bound cyanine dimers](#)

[Figure 8.2 — T3602; TO-PRO-1 iodide \(515/531\)](#)

[Figure 8.3 — T3600; TOTO-1 iodide \(514/533\)](#)

[Figure 8.4 — Lambda bacteriophage *Hind*III fragments. POPO™-1, BOBO™-1, YOYO\(R\)-1, TOTO\(R\)-1, JOJO™-1, POPO™-3, LOLO™-1, BOBO™-3, YOYO\(R\)-3 and TOTO\(R\)-3 nucleic acid stains.](#)

[Figure 8.5 — Relaxation of a single DNA molecule. YOYO\(R\)-1 iodide.](#)

[Figure 8.6 — NMR solution structure of the TOTO-1 dye bound to DNA](#)

[Figure 8.7 — Bovine pulmonary artery endothelial cells \(BPAEC\). MitoTracker\(R\) Red CMXRos, SYTOX\(R\) Green nucleic acid stain, biotin-XX goat anti-mouse IgG antibody and Cascade Blue\(R\) NeutrAvidin biotin-binding protein.](#)

[Figure 8.8 — Zebrafish gut. BODiPY\(R\) TR-X phalloidin and SYTOX\(R\) Green nucleic acid stain.](#)

[Figure 8.9 — *Micrococcus luteus* stained with SYTOX\(R\) Blue nucleic acid stain.](#)

[Figure 8.10 — SYTOX Orange/DNA](#)

[Figure 8.11 — Human neutrophil. SYTO\(R\) 13 live-cell nucleic acid stain.](#)

[Figure 8.12 — Relative absorbance and fluorescence emission spectra of SYTO RNASelect green-fluorescent cell stain in the presence of *Escherichia coli* RNA or *E. coli* DNA, or in buffer alone](#)

[Figure 8.13 — E1305; DISCONTINUED ethidium bromide](#)

[Figure 8.14 — P1304MP; propidium iodide](#)

[Figure 8.15 — *Drosophila* ovarian egg chamber. Propidium iodide.](#)

[Figure 8.16 — Fluorescence emission spectra of other DNA-bound dyes](#)

[Figure 8.17 — H7593; hexidium iodide](#)

[Figure 8.18 — E1169; ethidium homodimer-1 \(EthD-1\)](#)

[Figure 8.19 — E3599; ethidium homodimer-2 \(EthD-2\)](#)

[Figure 8.20 — Ethidium homodimer-1/DNA](#)

[Figure 8.21 — Fluorescence emission spectra of DNA-bound SYTOX Green stain and EthD-1](#)

[Figure 8.22 — E1374; ethidium monoazide bromide \(EMA\)](#)

[Figure 8.23 — A1301; acridine orange](#)

[Figure 8.24 — A666; acridine homodimer \(bis-\(6-chloro-2-methoxy-9-acridinyl\)spermine\)](#)

[Figure 8.25 — A1324; 9-amino-6-chloro-2-methoxyacridine \(ACMA\)](#)

[Figure 8.26 — H1398; Hoechst 33258, pentahydrate \(bis-benzimide\)](#)

[Figure 8.27 — H1399; Hoechst 33342, trihydrochloride, trihydrate](#)

[Figure 8.28 — Hoechst 33258/DNA](#)

[Figure 8.29 — D1306; 4',6-diamidino-2-phenylindole, dihydrochloride \(DAPI\)](#)

[Figure 8.30 — X-ray crystal structure of DAPI bound to DNA](#)

[Figure 8.31 — A1310; 7-aminoactinomycin D \(7-AAD\)](#)

[Figure 8.32 — 7-aminoactinomycin D/DNA](#)

[Figure 8.33 — Wheat root tips in seven stages of the cell cycle. 7-aminoactinomycin D.](#)

[Figure 8.34 — H22845; hydroxystilbamidine, methanesulfonate](#)

[Figure 8.35 — Fluorescence spectra of hydroxystilbamidine bound to different forms of DNA](#)

[Figure 8.36 — L7595; LDS 751](#)

[Figure 8.37 — Pyramidal cells of a mouse hippocampus and dentate gyrus. NeuroTrace\(R\) green fluorescent Nissl stain.](#)

[Figure 8.38 — Structure of ChromaTide BODIPY FL-14-dUTP](#)

[Figure 8.39 — Structure of ChromaTide Alexa Fluor 488-7-OBEA-dCTP](#)

[Figure 8.40 — F32767; fluorescein-aha-dUTP](#)

[Figure 8.41 — Expression of snail RNA in an early-stage fruit fly embryo visualized by FISH.](#)

[Figure 8.42 — Schematic diagram of the labeling method provided in our ULYSIS Nucleic Acid Labeling Kits](#)

[Figure 8.43 — Nucleic acid labeling method provided in our ULYSIS Nucleic Acid Labeling Kits](#)

[Figure 8.44 — Centromere probes to chromosomes 1, 15 and 17 labeled with ULYSIS Alexa Fluor 546, Alexa Fluor 594 and Oregon Green 488 Nucleic Acid Labeling Kits](#)

[Figure 8.45 — Schematic diagram of the labeling method provided in our ARES DNA Labeling Kits](#)

[Figure 8.46 — Human metaphase chromosome spreads.](#)

[Figure 8.47 — A21664; aminoallyl dUTP \(5-\(3-aminoallyl\)-2'-deoxyuridine 5'-triphosphate, trisodium salt\)](#)

[Figure 8.48 — QSY 7 dye-labeled oligonucleotide/pH 7.5](#)

[Figure 8.49 — Schematic diagram of labeling amine-modified oligodeoxynucleotides with succinimidyl ester dyes](#)

[Figure 8.50 — Fluorescence emission spectra of single-stranded DNA labeled with the Alexa Fluor 647 dye or Cy5 dye](#)

[Figure 8.51 — Fluorescence quenching of 5'-tetramethylrhodamine-labeled M13 primers by nonfluorescent dyes attached at the 3'-end](#)

[Figure 8.52 — DNA, RNA or protein quantitation with the Quant-iT Assay Kits](#)

[Figure 8.53 — DNA selectivity and sensitivity of the Quant-iT DNA High-Sensitivity assay](#)

[Figure 8.54 — DNA selectivity and sensitivity of the Quant-iT DNA Broad-Range assay](#)

[Figure 8.55 — RNA selectivity and sensitivity of the Quant-iT RNA assay](#)

[Figure 8.56 — Fluorescence enhancement of the PicoGreen reagent upon binding to dsDNA, ssDNA and RNA](#)

[Figure 8.57 — Quick and easy protocol of our nucleic acid quantitation assays](#)

[Figure 8.58 — Linear quantitation of calf thymus DNA using the PicoGreen dsDNA quantitation reagent](#)

[Figure 8.59 — Linear quantitation of DNA from 0 to 50 ng/mL using the PicoGreen dsDNA quantitation reagent](#)

[Figure 8.60 — A RediPlate 96 microplate](#)

[Figure 8.61 — Linear quantitation of a synthetic 24-mer using the OliGreen ssDNA quantitation reagent](#)

[Figure 8.62 — Linear quantitation of ribosomal RNA using the RiboGreen RNA quantitation reagent](#)

[Figure 8.63 — Dynamic range and sensitivity of the RediPlate 96 RiboGreen RNA quantitation assay](#)

[Figure 8.64 — Linear range and sensitivity of the RiboGreen RNA-specific quantitation assay](#)

[Figure 8.65 — Marine viruses, bacteria and a diatom. SYBR\(R\) Green I nucleic acid stain.](#)

[Figure 8.66 — Schematic representation of real-time PCR with the SYBR Green I dye](#)

[Figure 8.67 — Schematic representation of real-time PCR with quencher-labeled primers and the SYBR Green I dye](#)

[Figure 8.68 — DNA molecular weight ladders stained with SYBR\(R\) Green I nucleic acid gel stain.](#)

[Figure 8.69 — 1% Agarose gel containing 16S and 23S ribosomal RNA \(rRNA\). SYBR\(R\) Green II RNA gel stain.](#)

[Figure 8.70 — Comparison of glyoxalated RNA stained with ethidium bromide or SYBR Gold nucleic acid stain](#)

[Figure 8.71 — Comparison of sensitivity achieved using SYBR Gold stain and silver stain](#)

[Figure 8.72 — Fluorescence excitation and emission spectra of SYBR Gold nucleic acid gel stain bound to dsDNA](#)

[Figure 8.73 — Single-strand conformation polymorphism \(SSCP\) in exon 1 of human K-ras. SYBR\(R\) Gold nucleic acid gel stain.](#)

[Figure 8.74 — Comparison of dsDNA detection in native gels using SYBR Green I gel stain and ethidium bromide](#)

[Figure 8.75 — Agarose gel containing camptothecin-treated HL-60 cells. SYBR\(R\) Green I nucleic acid gel stain.](#)

[Figure 8.76 — Comparison of single-stranded oligonucleotide detection using SYBR Green I gel stain and ethidium bromide](#)

[Figure 8.77 — Fluorescence excitation and emission spectra of dsDNA-bound SYBR Green I gel stain](#)

[Figure 8.78 — Comparison of RNA detection in nondenaturing gels using SYBR Green II RNA gel stain and ethidium bromide](#)

[Figure 8.79 — Normalized fluorescence excitation and emission spectra of the SYBR Safe DNA gel stain](#)

[Figure 8.80 — DNA fragments stained with SYBR\(R\) Safe DNA gel stain.](#)

[Figure 8.81 — Human metaphase chromosomes hybridized to chromosome paints.](#)

[Figure 8.82 — Karyotype of human metaphase chromosomes using RxFISH.](#)

[Figure 8.83 — Visualization of chromosome 2 using a paint probe labeled with the ULYSIS Alexa Fluor\(R\) 546 Nucleic Acid Labeling Kit.](#)

[Figure 8.84 — Human metaphase chromosomes. ULYSIS Oregon Green\(R\) and Alexa Fluor\(R\) 594 Nucleic Acid Labeling Kits, and DAPI.](#)

[Figure 8.85 — Human metaphase chromosomes hybridized to centromere probes labeled with ChromaTide\(R\) OBEA-dCTP nucleotides.](#)

[Figure 8.86 — Fluorescence *in situ* hybridization \(FISH\) mapping of a BAC clone on human metaphase chromosomes](#)

[Figure 8.87 — Fluorescence *in situ* hybridization \(FISH\) to human metaphase chromosomes using Alexa Fluor conjugates for detection](#)

[Figure 8.88 — Schematic representation of *in situ* hybridization detection using amplification with biotin and a labeled streptavidin](#)

[Figure 8.89 — Schematic representation of *in situ* hybridization detection using amplification with an anti-dye antibody](#)

[Figure 8.90 — Amplification of FISH signals using the Alexa Fluor\(R\) 488 Signal-Amplification Kit for Fluorescein- and Oregon Green\(R\) Dye-Conjugated Probes.](#)

[Figure 8.91 — Schematic representation of mRNA *in situ* hybridization detection using tyramide signal amplification \(TSA\)](#)

[Figure 8.92 — Fluorescence *in situ* hybridization detected by tyramide signal amplification](#)

[Figure 8.93 — Schematic representation of mRNA *in situ* hybridization detection using the ELF technology](#)

[Figure 8.94 — Transformed mouse fibroblasts \(CRE BAG 2 cells\). ELF\(R\) 97 mRNA *In Situ* Hybridization Kit and DAPI](#)

[Figure 8.95 — HybriWell hybridization sealing system](#)

[Figure 8.96 — DNA microarray hybridized to DNA labeled using the ARES™ DNA Labeling Kits.](#)

[Figure 8.97 — Comparison of the absorption and fluorescence emission spectra of the Alexa Fluor 555 and Cy3 dyes](#)

[Figure 8.98 — Comparison of the fluorescence spectra of the Alexa Fluor 647 and Cy5 dyes](#)

[Figure 8.99 — ULYSIS reagent-labeled RNA hybridized to a microarray.](#)

[Figure 8.100 — Microarray hybridized to Alexa Fluor 546 dye-labeled cDNA](#)

[Figure 8.101 — Microarray expression analysis using an Alexa Fluor 546 dye-labeled dendrimer](#)

[Figure 8.102 — DNA microarrays stained with nucleic acid stains for quality control](#)

[Figure 8.103 — Panomer 9 oligodeoxynucleotides for quality control of microarray spotting](#)

[Figure 8.104 — Schematic diagram of method used in the PARAGON DNA Microarray QC Stain Kit](#)

[Figure 8.105 — Schematic diagram of the DNA microarray control slide from the PARAGON DNA Microarray QC Stain Kit and the PARAGON DNA Microarray QC Hybridization Kits](#)

[Figure 8.106 — Typical results obtained with the PARAGON Genomic DNA Hybridization Test Kit #2](#)

[Figure 8.107 — Typical results obtained with the PARAGON dye ratio calibration slide](#)

[Figure 8.108 — DNA stained with SYBR DX DNA blot stain following Southern transfer](#)

[Figure 8.109 — DNA microarray hybridized to Panomer 9 random oligodeoxynucleotides](#)

[Figure 8.110 — Microarray spotted with the Alexa Fluor 546 and Alexa Fluor 647 Panomer 9 random oligodeoxyribonucleotides](#)

[Figure 8.111 — Secure-Seal gasket for use with microarray hybridizations](#)

[Figure 8.112 — HybriSlip covers for hybridization](#)

[Figure 8.113 — Schematic representation of molecular beacons](#)

[Figure 8.114 — Schematic representation of wavelength-shifting molecular beacons](#)

[Figure 8.115 — Schematic representation of real-time PCR with TaqMan primers](#)

[Figure 8.116 — Schematic representation of real-time PCR with Scorpion primers](#)

[Figure 8.117 — Schematic representation of real-time PCR with UniPrimers](#)

[Figure 8.118 — Four-color image of mouse intestine cryosection.](#)

[Figure 8.119 — SYTOX Blue/DNA](#)

Figure 8.120 — Sexual reproduction of *Tetrahymena thermophila*. Texas Red(R)-X goat anti-mouse IgG (H+L) conjugate and SYTOX(R) Green nucleic acid stain.

Figure 8.121 — Macrophages treated with interleukin-4. Rhodamine phalloidin and YO-PRO(R)-1 iodide.

Figure 8.122 — Bovine pulmonary artery endothelial cells. SYTOX(R) Green nucleic acid stain and BODIPY(R) TR-X phalloidin.

Figure 8.123 — Human cheek epithelial cells. Alexa Fluor(R) 350 wheat germ agglutinin and SYTOX(R) Green nucleic acid stain.

Figure 8.124 — Nuclear deformation of an apoptotic cell visualized with the SYTOX(R) Green dye.

Figure 8.125 — Mouse brain section stained with NeuroTrace(R) 435/455 blue-fluorescent Nissl stain and counterstained with nuclear yellow.

Figure 8.126 — Binucleate bovine pulmonary artery endothelial cell. Anti- α -tubulin, biotin-XX conjugate; Alexa Fluor(R) 568 streptavidin and nuclear yellow.

Figure 8.127 — T1323; true blue chloride

Figure 8.128 — Zebrafish retina. TSA Kit #2. Alexa Fluor(R) 350 wheat germ agglutinin conjugate and TOTO(R)-3 nucleic acid stain.

Figure 8.129 — TOTO-3/DNA

Figure 8.130 — TO-PRO-3/DNA

Figure 8.131 — *Solemya reidi* gill filaments. The SYTOX(R) Green nucleic acid stain, propidium iodide or DAPI, which are all components of the FISH Counterstain Kit #2.

Figure 8.132 — Freshwater snail embryo. SYTO(R) Blue fluorescent nucleic acid stain, a component of the SYTO(R) Blue Fluorescent Stain Sampler Kit.

Figure 8.133 — Polymorphonuclear leukocytes. SYTO(R) 13 live-cell nucleic acid stain.

Figure 8.134 — Chromatin in PtK2 cells. DAPI.

Figure 8.135 — *Drosophila* embryos. Lissamine rhodamine B sulfonyl chloride and SYTOX(R) Green nucleic acid stain.

Figure 8.136 — Bovine pulmonary artery endothelial cells. Anti- α -tubulin, biotin-XX goat anti-mouse IgG antibody, Marina Blue(R) streptavidin and TO-PRO(R)-3.

Figure 8.137 — Mitotic spindles from sea urchin eggs. YOYO(R)-1 iodide and Texas Red(R) goat anti-rat IgG (H+L) antibody.

Figure 8.138 — Using the YOYO-1 dye to follow cell division in a sea urchin egg

Figure 8.139 — Centromeres in human metaphase chromosome spread. Fluorescein streptavidin, biotinylated anti-fluorescein and propidium iodide.

Figure 8.140 — Human metaphase chromosomes. DAPI.

Figure 8.141 — Human metaphase chromosomes. Chromosome Banding Kit #2 *SYTOX(R) Green nucleic acid stain/Methyl Green*.

Figure 8.142 — Mouse brain section stained with NeuroTrace(R) 515/535 yellow-fluorescent Nissl stain and counterstained with nuclear yellow.

Figure 8.143 — Mouse cerebellum section. NeuroTrace(R) 530/615 red-fluorescent Nissl stain and nuclear yellow.

Figure 8.144 — T2 phage genomic DNA. YOYO(R)-1 nucleic acid stain.

Figure 8.145 — Electrophoretic bandshift assays using SYBR Green I nucleic acid gel stain

Figure 8.146 — Titration of *lac* operator DNA with *lac* repressor protein.

Figure 8.147 — HL-60 cells in a comet assay. SYBR(R) Green I nucleic acid gel stain.

[Figure 8.148 — Aldehyde-reactive probe \(ARP\) used to detect DNA damage](#)

[Figure 8.149 — Schematic diagram of the mechanism used in the EnzChek Reverse Transcriptase Assay Kit](#)

[Figure 8.150 — Detection of HIV reverse transcription using the EnzChek Reverse Transcriptase Assay Kit](#)


List of Technical Notes and Product Highlights

[Note 8.1 — Product Highlight: Alexa Fluor Dyes for Labeling Nucleic Acids](#)

[Note 8.2 — Product Highlight: SYBR Green I and SYBR Green II Dyes for Quantitating Nucleic Acids on Plastic Wrap or Paraffin Sheets](#)

[Note 8.3 — Product Highlight: SYBR Safe DNA Gel Stain](#)

[Note 8.4 — Product Highlight: Procedure for Staining Microarrays with SYBR Green II Dye](#)

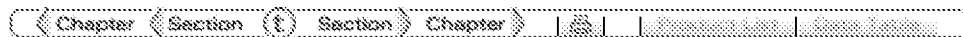

[Site Map](#) | [Change Country](#) | [Contact Us](#) | [Privacy Policy](#) | [Terms & Conditions](#) | [Online Price Policy](#) | [Feedback](#) | [Careers](#)

Copyright © 2008 Invitrogen Corporation. All rights reserved.