



[My Account](#) | [Order Status](#) | [Quick Order Worldwide](#)



- [Products & Services](#)
  - [Applications](#)
    - [Nucleic Acid Purification & Analysis](#)
    - [Nucleic Acid Amplification & Expression Profiling](#)
    - [Cloning](#)
    - [Proteins, Expression, Isolation & Analysis](#)
    - [RNAi, Epigenetics & Gene Regulation](#)
    - [Cell Culture](#)
    - [Cell & Tissue Analysis](#)
    - [Drug Discovery](#)
    - [Clinical & Diagnostic Applications](#)
    - [Bioproduction](#)
    -
  - [Product Types](#)
    - [Antibodies & Secondary Detection Reagents](#)
    - [Beads & Microspheres](#)
    - [Buffers & Chemicals](#)
    - [Cell Culture & Transfection Reagents](#)
    - [Columns, Resins & Plates](#)
    - [Competent Cells & Strains](#)
    - [Controls, Standards & Ladders](#)
    - [DNA Vectors, Clones & Libraries](#)
    - [Eukaryotic Cells](#)
    - [Gels, Fractionation Strips, & Membranes](#)
    - [Instruments, Software, & Accessories](#)
    - [Kits & Assays](#)
    - [Labeling & Detection Products](#)
    - [Literature](#)
    - [Primers, Oligos, & Nucleotides](#)
    - [Proteins, Enzymes, & Peptides](#)
    - [Slide & Chip Arrays](#)
    -
  - [Services](#)
    - [Discovery Research](#)
    - [Production & Process Development](#)
    - [OEM Services](#)
    - [Cell Culture Services™](#)
    - [Licenses, Terms, and Conditions](#)
    -
  - [Promotions](#)
  - [New Products](#)
  - [Sample Offers](#)
  - [Check Price Tool](#)
  - [Quick Order Tool](#)
  - [Supply Center](#)
    - [Request Information](#)

- [SCMS Account Registration](#)
    - [SCMS Admin Login](#)
    -
  - [eProcurement](#)
  -
- 
- [Applications](#)
  - [Nucleic Acid Purification & Analysis](#)
    - [DNA Purification](#)
    - [RNA Purification](#)
    - [Viral RNA/DNA Purification](#)
    - [Nucleic Acid Quantitation](#)
    - [Nucleic Acid Gel Electrophoresis](#)
    - [Automated Nucleic Acid Purification](#)
    -
  - [Nucleic Acid Amplification & Expression Profiling](#)
    - [Oligonucleotide Design](#)
    - [Real-Time PCR](#)
    - [PCR](#)
    - [Genomic Profiling - Array Labeling & Genotyping](#)
    - [Reverse Transcription & cDNA Synthesis](#)
    - [Gene Expression Profiling - Array Labeling & SAGE](#)
    - [RT-PCR](#)
    - [DNA Sequencing](#)
    - [Nucleic Acid Labeling](#)
    -
  - [Cloning](#)
    - [PCR Cloning](#)
    - [PCR](#)
    - [Gateway® Cloning](#)
    - [Transformation](#)
    - [Restriction & Cloning Enzymes](#)
    - [Oligonucleotide Design](#)
    - [cDNA Libraries & Library Construction](#)
    - [Cloning Services](#)
    - [Mutagenesis](#)
    - [Premade Clones](#)
    - [Microbiological Culture](#)
    - [Vector NTI® Software](#)
    -
  - [Proteins, Expression, Isolation & Analysis](#)
    - [Antibodies and Purified Proteins](#)
    - [ELISA](#)
    - [Luminex® Assays](#)
    - [Protein Expression Systems & Vectors](#)
    - [Protein Quantitation](#)
    - [Protein-Protein Interactions](#)
    - [Protein Sample Preparation & Protein Isolation](#)
    - [Mass Spectrometry](#)
    - [Enzyme & Protein Activity Assays](#)
    - [Transfection & Selection](#)
    - [Protein Gel Electrophoresis](#)
    - [Western Blotting](#)

- Biomarker Discovery
- 
- RNAi, Epigenetics & Gene Regulation
  - RNAi
  - miRNA Profiling
  - Methylation Analysis
  - Transcription Factor Binding Assays
  - Chromatin Immunoprecipitation (ChIP)
  -
- Cell Culture
  - Mammalian Cell Culture
  - Stem Cell Research
  - Transfection & Selection
  - Insect Cell Culture
  - 3D Cell Culture
  - Cell Counting, Viability & Cryopreservation
  - Cell Isolation & Expansion
  - Primary Cell Culture
  -
- Cell & Tissue Analysis
  - Antibodies & Secondary Detection
  - Cell Tracing, Tracking
  - Cell Isolation & Expansion
  - Cellular Imaging
  - Stem Cell Differentiation
  - Flow Cytometry
  - Cell Viability, Proliferation & Function
  - Immunoassays
  - Cell Structure
  - Labeling Chemistry
  - Quantum Dots & Microspheres
  - Fluorescent Microspheres
  - Microbiological Analysis
  -
- Drug Discovery
  - Target & Lead Identification & Validation
  - Drug Metabolism & Safety
  -
- Clinical & Diagnostic Applications
  - Transplant Diagnostics
  - Clinical Cell Isolation & Expansion
  - Anatomical Pathology
  - Cytogenetics
  - BioDefense & Pathogen Detection
  - Environmental Diagnostics
  - Bead-Based IVD Assays
  -
- BioProduction
  - Protein & Antibody Production
  - Virus & Vaccine Production
  - Custom Media Products
  - Bioproduction Services
  - Basic Cell Culture

- 
- Brands
  - Dynal®
    - Dynabeads® Products & Technology
    - The History of Dynal® and Biomagnetic Separation
    - Bead Types
    - Magnets
  - Gibco
    - New Technologies
    - General Information
    - Reference Guide
    - Press Releases
  - Molecular Probes®
    - The Handbook
    - Literature
    - Subscribe to ProbesOnline Newsletter
    - Legal Notices and Trademark Attributions
    - Licensing and Business Development
  - Product Brands
  -
- 
- Guides
  - LINNEA™ Communities
    - Vector NTI® Community
  - LINNEA™ Designers
    - BLOCK-iT™ RNAi Designer
    - D-Lux™ Designer
    - OligoPerfect™ Designer
    - PeptideSelect™
    - VectorDesigner™
  - LINNEA™ Genes
  - LINNEA™ Guide to Antibodies
  - LINNEA™ Guide to Clones
    - Ultimate™ ORF Clones
    - Clone Search
  - LINNEA™ Guide to Primers
    - D-Lux™ Assay Search
    - D-Lux™ Designer
    - OligoPerfect™ Designer
  - LINNEA™ Guide to RNAi
    - BLOCK-iT™ RNAi Express
    - BLOCK-iT™ RNAi Designer
  - LINNEA™ Pathways

- [LINNEA™ Protocols](#)
  -
- 
- [Support](#)
  - [Order Support - FAQ](#)
  - [Product Support - FAQ](#)
  - [Scientific Tools & Tips](#)
  - [Newsletters & Journals](#)
    - [Cell Culture Archive](#)
    - [Benchtopics Archive](#)
    - [Focus Journal Archive](#)
    - [Expressions Archive](#)
    - [Postings Archive](#)
    - [Quest Archive](#)
    -
  - [Research Tools](#)
    - [KD calculator](#)
    - [Image Gallery](#)
    - [Base-Dye Ratio Calculator](#)
    - [Fluorescence SpectraViewer](#)
    -
  - [Online Seminars](#)
  - [Product Technical Resources](#)
    - [Product Spectra](#)
    - [Product Structures](#)
    - [Citations and Bibliographic References Search](#)
    - [Cell Lines Advanced Search](#)
    - [Cell Lines Browse](#)
    -
  - [Tutorials](#)
  - [Search Tips](#)
  -
- 
- [Corporate](#)
  - [Who We Are](#)
    - [Our Quest](#)
    - [Company Fact Sheet](#)
    - [Areas of Focus](#)
    - [Global Locations & Directions](#)
    - [Awards & Recognition](#)
    - [Customer Testimonials](#)
    -
  - [Corporate Citizenship](#)
    - [Community Involvement](#)
    - [Invitrogen Foundation](#)
    - [Environmental Responsibility](#)
    -
  - [Careers](#)
    - [Job Listings](#)
    - [Inside Invitrogen](#)
    - [Meet Our Employees](#)
    - [Benefits](#)
    -

- [Media Center](#)
  - [Press Releases](#)
  - [Image Gallery](#)
  - [Videos & Podcasts](#)
  - [Media Contacts](#)
  - [News Subscription](#)
  -
- [Investor Relations](#)
  - [All SEC Filings](#)
  - [Calendar of Events](#)
  - [Corporate Governance](#)
  - [Investor FAQs](#)
  - [Financials](#)
  - [Investor Kits](#)
  - [Presentations](#)
  - [Press Releases](#)
  - [Shareholder Information](#)
  - [Stock Information](#)
  - [Web Cast Archive](#)
  -
- 

**Search**

All Results

[Search](#)

[View Cart](#)



0 Items

[LOGIN](#)

## Molecular Probes

- [The Handbook](#)
  - [Recent Updates](#)
  - [Technical Notes and Product Highlights](#)
  - [List of Tables](#)
  - [List of Figures](#)
- [Literature](#)
- [Subscribe to ProbesOnline Newsletter](#)
- [Subscribe to Our Mailing List](#)
- [Legal Notices and Trademark Attributions](#)
- [Licensing and Business Development](#)

Have a question or  
just need help?

Contact us for  
Service and Support




[Home](#) > [Brands](#) > [Molecular Probes](#) > [Literature](#) > [The Handbook](#)

Chapter Section Section Chapter

## Chapter 8 — Nucleic Acid Detection and Genomics Technology

### Section 8.1 — Nucleic Acid Stains

 [Printer-Friendly Format](#)

- [Properties of Cyanine Dyes](#)
- [Premier Cyanine Dyes for Ultrasensitive Nucleic Acid Detection and Quantitation](#)
- [Cell-Impermeant Cyanine Dimers: The TOTO Family of Dyes](#)
  - [High Affinity for Nucleic Acids](#)
  - [High Fluorescence Enhancements and High Quantum Yields upon Binding to Nucleic Acids](#)
  - [Modifying the Dimers Creates Compounds with Different Spectral Characteristics](#)
  - [Binding Modes of the Cyanine Dimers](#)
  - [Working with Cyanine Dimers](#)
- [Cell-Impermeant Cyanine Monomers: The TO-PRO Family of Dyes](#)
  - [Spectral Characteristics of the Cyanine Dye Monomers](#)
  - [Working with Cyanine Monomers](#)
- [Cell-Impermeant SYTOX Dyes for Dead-Cell Staining](#)
  - [SYTOX Green Stain](#)
  - [SYTOX Blue Stain](#)
  - [SYTOX Orange Stain](#)
- [Cell-Permeant Cyanine Dyes: The SYTO Nucleic Acid Stains](#)
  - [SYTO Nucleic Acid Stains for DNA and RNA](#)
  - [SYTO RNASelect Green-Fluorescent Cell Stain](#)
- [Chemically Reactive Cyanine Dyes](#)
- [RediPlate 96 Nucleic Acid Stain Sampler Microplate](#)
- [Phenanthridines and Acridines: Classic Intercalating Dyes](#)
  - [Cell-Impermeant Ethidium Bromide and Propidium Iodide](#)
  - [Cell-Permeant Hexidium Iodide](#)
  - [Cell-Permeant Dihydroethidium \(Hydroethidine\)](#)
  - [High-Affinity Ethidium Homodimers](#)
  - [Ethidium Monoazide: A Photocrosslinking Reagent](#)
  - [Acridine Orange: A Dual-Fluorescence Nucleic Acid Stain](#)
  - [AT-Selective Acridine Homodimer](#)
  - [AT-Selective ACMA](#)
- [Indoles and Imidazoles: Classic Minor Groove-Binding Dyes](#)
  - [DNA-Selective Hoechst Dyes](#)
  - [AT-Selective DAPI](#)
- [Other Nucleic Acid Stains](#)
  - [7-Aminoactinomycin D and Actinomycin D: Fluorescent Intercalators](#)
  - [Multicolor Hydroxystilbamidine](#)
  - [Long-Wavelength LDS 751](#)
  - [NeuroTrace Fluorescent Nissl Stains](#)
- [Data Table](#)
- [Product List](#)

### Section 8.2 — Labeling Oligonucleotides and Nucleic Acids

 [Printer-Friendly Format](#)

- [ChromaTide Nucleotides](#)
  - [Structures of the ChromaTide Nucleotides](#)
  - [Fluorescent ChromaTide Nucleotides](#)
  - [ChromaTide Dinitrophenyl \(DNP\)-11-dUTP](#)
  - [Using ChromaTide Nucleotides in Enzymatic Labeling Methods](#)
- [Amine-Modified Nucleotides](#)
  - [Unlabeled and Labeled aha-dUTP and aha-dCTP](#)
  - [Aminoallyl UTP and Aminoallyl dUTP](#)
  - [Alexa Fluor Amine-Reactive Dye Decapacks for Labeling Amine-Modified DNA and RNA](#)


- 5-Bromo-2'-Deoxyuridine, 5-Bromo-dUTP (BrdUTP) and 5-Bromo-UTP (BrUTP)
- ULYSIS Nucleic Acid Labeling Kits
- ARES DNA Labeling Kits
- Labeled Oligonucleotides
- Labeling Amine- and Thiol-Modified Oligonucleotides
  - Alexa Fluor Oligonucleotide Amine Labeling Kits
  - Alexa Fluor Amine-Reactive Dye Decapacks
  - Other Reactive Fluorescent Dyes and Quenchers
  - Dyes for Sequencing Applications
- Labeling Phosphate-Modified Oligonucleotides
- Other Chemical Labeling Methods for Nucleic Acids
  - Labeling Cytidine Residues
  - Labeling the 3'-Terminus of RNA
  - Labeling Abasic Sites with ARP
  - Specialized Methods for Nucleic Acid Modification
- ElutaTube Microdialysis Vials
- [Product List](#)

## Section 8.3 — Nucleic Acid Detection and Quantitation in Solution



- Nucleic Acid Detection and Quantitation in Solution
- QuantiT Assay Kits for DNA and RNA
  - QuantiT DNA Assay Kits
  - QuantiT RNA Assay Kit
- PicoGreen dsDNA Quantitation Assay
  - PicoGreen dsDNA Quantitation Reagent and Kits
  - RediPlate 96 PicoGreen dsDNA Quantitation Kits and Microplate
- OliGreen ssDNA Quantitation Reagent and Kit
- RiboGreen RNA Quantitation Assay
  - RiboGreen RNA Quantitation Reagent and Kit
  - RediPlate 96 RiboGreen RNA Quantitation Kit
  - RiboGreen RNA-Specific Quantitation Kit with DNase I
- Other Stains for Nucleic Acid Quantitation in Solution
  - Cyanine Dyes and Phenanthridine Dyes for Nucleic Acid Quantitation in Solution
  - Hoechst 33258 Dye for Quantitating DNA in Solution
- Real-Time Quantitative PCR Using the SYBR Green I Nucleic Acid Gel Stain
- A Useful Technique to Reduce Background Fluorescence in Nucleic Acid Detection
- [Product List](#)

## Section 8.4 — Nucleic Acid Detection and Quantitation in Electrophoretic Gels and Capillaries

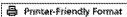

- Nucleic Acid Detection in Gels
  - SYBR Gold Nucleic Acid Gel Stain
  - SYBR Green I Nucleic Acid Gel Stain
  - SYBR Green II RNA Gel Stain
  - SYBR Green Nucleic Acid Gel Stains: Special Packaging and a Starter Kit
  - SYBR Safe DNA Gel Stain
  - Ethidium Bromide
  - Cyanine Monomers for Staining DNA in Electrophoretic Gels
  - Cyanine and Ethidium Dimers for Staining DNA Prior to Electrophoresis
  - Electrophoretic Mobility-Shift (Bandshift) Assay (EMSA) Kit
  - Other Nucleic Acid Stains for Gel-Staining Applications
- Safe Imager Blue-Light Transilluminator
- Capillary Electrophoresis and Channel Electrophoresis
  - Capillary Electrophoresis
- [Product List](#)

## Section 8.5 — Detecting Nucleic Acid Hybridization




- Principles of Fluorescence *In Situ* Hybridization (FISH)
- Probe Preparation for FISH



- Enzymatic Incorporation of ChromaTide Labeled Nucleotides
  - Aminoallyl dUTP Labeling Using the ARES Kits
  - Direct Chemical Labeling of Nucleic Acids Using the ULYSIS Kits
  - Labeling Oligonucleotides
- Signal Amplification for FISH Using Secondary Detection Reagents
  - Biotin and Streptavidin
  - Signal Amplification Using Anti-Dye Antibodies
  - Tyramide Signal Amplification (TSA)
  - Signal Amplification Using Enzyme-Labeled Fluorescence (ELF) Technology
  - Signal Amplification Using Two Rounds of Amplification
  - Colorimetric Signal Amplification
- Counterstaining Chromosomes for FISH
- Gaskets for *In Situ* Hybridization Experiments
- Nucleic Acid Hybridization on Blots and Microarrays
- Labeling Nucleic Acid Samples for Microarray Experiments
  - Consistent, Uniform Sample Labeling with ARES DNA Labeling Kits
  - Alexa Fluor Reactive Dye Decapacks for Microarray Applications
  - Fast and Easy Direct Labeling with ULYSIS Nucleic Acid Labeling Kits
  - Fluorophore-Labeled Nucleotides
  - Labeled dT<sub>18</sub> and Random Oligodeoxynucleotides
  - Using Labeled Oligonucleotides as Biosensors on Solid Supports
- Secondary Detection for Signal Amplification
  - Detecting Biotin Labels with Fluorescent Streptavidin Conjugates
  - Signal Amplification with the ELF 97 Phosphatase Substrate
  - Chemiluminescent Detection on Blots
  - Other Types of Signal Amplification
- Tools and Techniques for Quality Control Testing and Normalization on Microarrays and Blots
  - PARAGON DNA Microarray QC Stain Kit
  - PARAGON DNA Microarray QC Hybridization Kits
  - PARAGON Genomic DNA Hybridization Test Kits
  - PARAGON Microarray Scanner Calibration Slide
  - PARAGON Dye Ratio Calibration Slide
  - Other Nucleic Acid Stains for Standardizing Microarrays
  - SYBR DX DNA Biot Stain
  - Panomer Random-Sequence Oligonucleotides
- Sample Quantitation
- Hybridization Chambers for Microarrays
- Solution-Based Hybridization Assays
  - Real-Time Quantitative PCR Using FRET and Quenching Techniques
  - Other Solution-Based Hybridization Assays Using FRET and Quenching Techniques
  - Other Solution-Based Hybridization Assays Using Fluorescent Nucleic Acid Stains
- [Product List](#)

## Section 8.6 — Nuclear and Chromosome Counterstaining and Nissl Stains

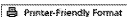

- Nuclear Counterstaining of Fixed Cells and Tissues
  - Blue-Fluorescent Counterstains
  - Green-Fluorescent Counterstains
  - Yellow-Fluorescent Counterstain
  - Orange- and Red-Fluorescent Counterstains
  - Long-Wavelength Nuclear Counterstains
  - SelectFX Nuclear Labeling Kit
  - Cytological Nuclear Counterstain Kit
- Nuclear Counterstains for Live Cells and Unfixed Tissues
  - Cell-Permeant Blue-Fluorescent Counterstains
  - Cell-Permeant Green-Fluorescent Counterstains
  - Cell-Permeant Orange- and Red-Fluorescent Counterstains
- Tracking Chromosomes through Mitosis
- Chromosome Counterstaining
  - Red-Fluorescent Counterstaining for Green-Fluorescent Probes
  - The Best Green-Fluorescent Counterstains
  - Blue-Fluorescent Counterstaining for Multicolor Labeling
  - Very Long-Wavelength Counterstains for Multicolor Labeling
- Chromosome Banding
  - SYTOX Green Nucleic Acid Stain
  - Acridine Homodimer
  - Other Dyes and Chromosome Banding Reagents
- NeuroTrace Fluorescent Nissl Stains
- [Product List](#)

## Section 8.7 — Analysis of DNA Structure, DNA Binding and DNA Damage

 [Printer-Friendly Format](#)

- Nucleic Acid Conformational Analysis
- Examining the Behavior of Single Nucleic Acid Molecules
- DNA Binding Assays
  - Electrophoretic Mobility-Shift (Bandshift) Assays
  - DNA Binding Assays in Solution
  - Selective Cleavage of Nucleic Acids with a Chemical Nuclease
- Assessing DNA Damage
  - Comet (Single-Cell Gel Electrophoresis) Assay to Detect Damaged DNA
  - TUNEL Assay for *In Situ* Detection of Fragmented DNA
  - Microplate Assays for DNA Damage
- Assays for Enzymes that Modify Nucleic Acids
  - Gel Assays for DNase Detection
  - Solution Assays for Nuclease Detection
  - An Assay for Reverse Transcriptase Activity
  - Telomerase
- [Product List](#)

## List of Tables

[Table 8.1 — Specialty nucleic acid reagents for molecular biology](#)

[Table 8.2 — Cell membrane-impermeant cyanine nucleic acid stains](#)

[Table 8.3 — Cell-permeant cyanine nucleic acid stains](#)

[Table 8.4 — Properties of classic nucleic acid stains](#)

[Table 8.5 — General characteristics of the dyes provided in the RediPlate 96 nucleic acid stain sampler microplate \(R32715\)](#)

[Table 8.6 — Characteristics of ChromaTide UTP nucleotides](#)

[Table 8.7 — Characteristics of ChromaTide dUTP, ChromaTide OBFA-dCTP, aha-dUTP and aha-dCTP labeled nucleotides](#)

[Table 8.8 — Spectral characteristics of the fluorescent dyes available in Molecular Probes' ULYSIS Nucleic Acid Labeling Kits](#)

[Table 8.9 — Spectral characteristics of the fluorescent dyes available in Molecular Probes' ARES DNA Labeling Kits](#)

[Table 8.10 — Oligonucleotide Amine Labeling Kits](#)

[Table 8.11 — Amine-reactive dyes for nucleic acid sequencing](#)

[Table 8.12 — Selection Guide for the Quant-iT Assay Kits](#)

[Table 8.13 — Effects of Contaminants on the PicoGreen Assay](#)

[Table 8.14 — Effects of Contaminants on the OliGreen Assay](#)

[Table 8.15 — Comparison of the RiboGreen RNA quantitation reagent and ethidium bromide](#)

[Table 8.16 — Effects of Contaminants on the RiboGreen Assay](#)

[Table 8.17 — Tools for hybridization experiments](#)

[Table 8.18 — Spectral characteristics of labeled oligonucleotides](#)

[Table 8.19 — Nucleic acid stains used for microarray quality control](#)

## List of Figures

[Figure 8.1 — Fluorescence emission spectra of DNA-bound cyanine dimers](#)

[Figure 8.2 — T3602; TO-PRO-1 iodide \(515/531\)](#)

[Figure 8.3 — T3600; TOTO-1 iodide \(514/533\)](#)

[Figure 8.4 — Lambda bacteriophage HindIII fragments. POPO™-1, BOBO™-1, YOYO\(R\)-1, TOTO\(R\)-1, JOJO™-1, POPO™-3, LOLO™-1, BOBO™-3, YOYO\(R\)-3 and TOTO\(R\)-3 nucleic acid stains.](#)

[Figure 8.5 — Relaxation of a single DNA molecule. YOYO\(R\)-1 iodide.](#)

[Figure 8.6 — NMR solution structure of the TOTO-1 dye bound to DNA](#)

[Figure 8.7 — Bovine pulmonary artery endothelial cells \(BPAEC\). MitoTracker\(R\) Red CMXRos, SYTOX\(R\) Green nucleic acid stain, biotin-XX goat anti-mouse IgG antibody and Cascade Blue\(R\) NeutrAvidin biotin-binding protein.](#)

[Figure 8.8 — Zebrafish gut. BODIPY\(R\) TR-X phalloidin and SYTOX\(R\) Green nucleic acid stain.](#)

[Figure 8.9 — \*Micrococcus luteus\* stained with SYTOX\(R\) Blue nucleic acid stain.](#)

[Figure 8.10 — SYTOX Orange/DNA](#)

[Figure 8.11 — Human neutrophil. SYTO\(R\) 13 live-cell nucleic acid stain.](#)

[Figure 8.12 — Relative absorbance and fluorescence emission spectra of SYTO RNaselect green-fluorescent cell stain in the presence of \*Escherichia coli\* RNA or \*E. coli\* DNA, or in buffer alone](#)

[Figure 8.13 — E1305; DISCONTINUED ethidium bromide](#)

[Figure 8.14 — P1304MP; propidium iodide](#)

[Figure 8.15 — \*Drosophila\* ovarian egg chamber. Propidium iodide.](#)

[Figure 8.16 — Fluorescence emission spectra of other DNA-bound dyes](#)

[Figure 8.17 — H7593; hexidium iodide](#)

[Figure 8.18 — E1169; ethidium homodimer-1 \(EthD-1\)](#)

[Figure 8.19 — E3599; ethidium homodimer-2 \(EthD-2\)](#)

[Figure 8.20 — Ethidium homodimer-1/DNA](#)

[Figure 8.21 — Fluorescence emission spectra of DNA-bound SYTOX Green stain and EthD-1](#)

[Figure 8.22 — E1374; ethidium monoazide bromide \(EMA\)](#)

[Figure 8.23 — A1301; acridine orange](#)

[Figure 8.24 — A666; acridine homodimer \(bis-\(6-chloro-2-methoxy-9-acridinyl\) spermine\)](#)

[Figure 8.25 — A1324, 9-amino-6-chloro-2-methoxyacridine \(ACMA\)](#)

[Figure 8.26 — H1398; Hoechst 33258, pentahydrate \(bis-benzimide\)](#)

[Figure 8.27 — H1399; Hoechst 33342, trihydrochloride, trihydrate](#)

[Figure 8.28 — Hoechst 33258/DNA](#)

[Figure 8.29 — D1306; 4',6-diamidino-2-phenylindole, dihydrochloride \(DAPI\)](#)

Figure 8.30 — X-ray crystal structure of DAPI bound to DNA

Figure 8.31 — A1310; 7-aminoactinomycin D (7-AAD)

Figure 8.32 — 7-aminoactinomycin D/DNA

Figure 8.33 — Wheat root tips in seven stages of the cell cycle, 7-aminoactinomycin D.

Figure 8.34 — H22845; hydroxystilbamidine, methanesulfonate

Figure 8.35 — Fluorescence spectra of hydroxystilbamidine bound to different forms of DNA

Figure 8.36 — L7595; LDS 751

Figure 8.37 — Pyramidal cells of a mouse hippocampus and dentate gyrus. NeuroTrace(R) green fluorescent Nissl stain.

Figure 8.38 — Structure of ChromaTide BODIPY FL-14-dUTP

Figure 8.39 — Structure of ChromaTide Alexa Fluor 488-7-OBEA-dCTP

Figure 8.40 — F32767; fluorescein-aha-dUTP

Figure 8.41 — Expression of snail RNA in an early-stage fruit fly embryo visualized by FISH.

Figure 8.42 — Schematic diagram of the labeling method provided in our ULYSIS Nucleic Acid Labeling Kits

Figure 8.43 — Nucleic acid labeling method provided in our ULYSIS Nucleic Acid Labeling Kits

Figure 8.44 — Centromere probes to chromosomes 1, 15 and 17 labeled with ULYSIS Alexa Fluor 546, Alexa Fluor 594 and Oregon Green 488 Nucleic Acid Labeling Kits

Figure 8.45 — Schematic diagram of the labeling method provided in our ARES DNA Labeling Kits

Figure 8.46 — Human metaphase chromosome spreads.

Figure 8.47 — A21864; aminoallyl dUTP (5-(3-aminoallyl)-2'-deoxyuridine 5'-triphosphate, trisodium salt)

Figure 8.48 — QSY 7 dye-labeled oligonucleotide/pH 7.5

Figure 8.49 — Schematic diagram of labeling amine-modified oligodeoxynucleotides with succinimidyl ester dyes

Figure 8.50 — Fluorescence emission spectra of single-stranded DNA labeled with the Alexa Fluor 647 dye or Cy5 dye

Figure 8.51 — Fluorescence quenching of 5'-tetramethylrhodamine-labeled M13 primers by nonfluorescent dyes attached at the 3'-end

Figure 8.52 — DNA, RNA or protein quantitation with the Quant-iT Assay Kits

Figure 8.53 — DNA selectivity and sensitivity of the Quant-iT DNA High-Sensitivity assay

Figure 8.54 — DNA selectivity and sensitivity of the Quant-iT DNA Broad-Range assay

Figure 8.55 — RNA selectivity and sensitivity of the Quant-iT RNA assay

Figure 8.56 — Fluorescence enhancement of the PicoGreen reagent upon binding to dsDNA, ssDNA and RNA

Figure 8.57 — Quick and easy protocol of our nucleic acid quantitation assays

Figure 8.58 — Linear quantitation of calf thymus DNA using the PicoGreen dsDNA quantitation reagent

Figure 8.59 — Linear quantitation of DNA from 0 to 50 ng/mL using the PicoGreen dsDNA quantitation reagent

Figure 8.60 — A RediPlate 96 microplate

Figure 8.61 — Linear quantitation of a synthetic 24-mer using the OliGreen ssDNA quantitation reagent

Figure 8.62 — Linear quantitation of ribosomal RNA using the RiboGreen RNA quantitation reagent

Figure 8.63 — Dynamic range and sensitivity of the RediPlate 96 RiboGreen RNA quantitation assay

Figure 8.64 — Linear range and sensitivity of the RiboGreen RNA-specific quantitation assay

Figure 8.65 — Marine viruses, bacteria and a diatom, SYBR(R) Green I nucleic acid stain.

Figure 8.66 — Schematic representation of real-time PCR with the SYBR Green I dye

Figure 8.67 — Schematic representation of real-time PCR with quencher-labeled primers and the SYBR Green I dye

Figure 8.68 — DNA molecular weight ladders stained with SYBR(R) Green I nucleic acid gel stain.

Figure 8.69 — 1% Agarose gel containing 16S and 23S ribosomal RNA (rRNA). SYBR(R) Green II RNA gel stain.

Figure 8.70 — Comparison of glyoxylated RNA stained with ethidium bromide or SYBR Gold nucleic acid stain

Figure 8.71 — Comparison of sensitivity achieved using SYBR Gold stain and silver stain

Figure 8.72 — Fluorescence excitation and emission spectra of SYBR Gold nucleic acid gel stain bound to dsDNA

Figure 8.73 — Single-strand conformation polymorphism (SSCP) in exon 1 of human K-ras. SYBR(R) Gold nucleic acid gel stain.

Figure 8.74 — Comparison of dsDNA detection in native gels using SYBR Green I gel stain and ethidium bromide

Figure 8.75 — Agarose gel containing camptothecin-treated HL-60 cells. SYBR(R) Green I nucleic acid gel stain.

Figure 8.76 — Comparison of single-stranded oligonucleotide detection using SYBR Green I gel stain and ethidium bromide

Figure 8.77 — Fluorescence excitation and emission spectra of dsDNA-bound SYBR Green I gel stain

Figure 8.78 — Comparison of RNA detection in nondenaturing gels using SYBR Green II RNA gel stain and ethidium bromide

Figure 8.79 — Normalized fluorescence excitation and emission spectra of the SYBR Safe DNA gel stain

Figure 8.80 — DNA fragments stained with SYBR(R) Safe DNA gel stain.

Figure 8.81 — Human metaphase chromosomes hybridized to chromosome paints.

Figure 8.82 — Karyotype of human metaphase chromosomes using RxFISH.

Figure 8.83 — Visualization of chromosome 2 using a paint probe labeled with the ULYSIS Alexa Fluor(R) 546 Nucleic Acid Labeling Kit.

Figure 8.84 — Human metaphase chromosomes. ULYSIS Oregon Green(R) and Alexa Fluor(R) 594 Nucleic Acid Labeling Kits, and DAPI.

Figure 8.85 — Human metaphase chromosomes hybridized to centromere probes labeled with ChromaTide(R) QBEA-dCTP nucleotides.

Figure 8.86 — Fluorescence *in situ* hybridization (FISH) mapping of a BAC clone on human metaphase chromosomes

Figure 8.87 — Fluorescence *in situ* hybridization (FISH) to human metaphase chromosomes using Alexa Fluor conjugates for detection

Figure 8.88 — Schematic representation of *in situ* hybridization detection using amplification with biotin and a labeled streptavidin

Figure 8.89 — Schematic representation of *in situ* hybridization detection using amplification with an anti-dye antibody

Figure 8.90 — Amplification of FISH signals using the Alexa Fluor(R) 488 Signal-Amplification Kit for Fluorescein- and Oregon Green(R) Dye-Conjugated Probes.

Figure 8.91 — Schematic representation of mRNA *in situ* hybridization detection using tyramide signal amplification (TSA)

Figure 8.92 — Fluorescence *in situ* hybridization detected by tyramide signal amplification

Figure 8.93 — Schematic representation of mRNA *in situ* hybridization detection using the ELF technology

Figure 8.94 — Transformed mouse fibroblasts (CRE BAG 2 cells). ELF(R) 97 mRNA *In Situ* Hybridization Kit and DAPI

Figure 8.95 — HybriWell hybridization sealing system

Figure 8.96 — DNA microarray hybridized to DNA labeled using the ARES™ DNA Labeling Kits.

Figure 8.97 — Comparison of the absorption and fluorescence emission spectra of the Alexa Fluor 555 and Cy3 dyes

Figure 8.98 — Comparison of the fluorescence spectra of the Alexa Fluor 647 and Cy5 dyes

Figure 8.99 — ULYSIS reagent-labeled RNA hybridized to a microarray.

Figure 8.100 — Microarray hybridized to Alexa Fluor 546 dye-labeled cDNA

Figure 8.101 — Microarray expression analysis using an Alexa Fluor 546 dye-labeled dendrimer

Figure 8.102 — DNA microarrays stained with nucleic acid stains for quality control

Figure 8.103 — Panomer 9 oligodeoxynucleotides for quality control of microarray spotting

Figure 8.104 — Schematic diagram of method used in the PARAGON DNA Microarray QC Stain Kit

Figure 8.105 — Schematic diagram of the DNA microarray control slide from the PARAGON DNA Microarray QC Stain Kit and the PARAGON DNA Microarray QC Hybridization Kits

Figure 8.106 — Typical results obtained with the PARAGON Genomic DNA Hybridization Test Kit #2

Figure 8.107 — Typical results obtained with the PARAGON dye ratio calibration slide

Figure 8.108 — DNA stained with SYBR DX DNA blot stain following Southern transfer

Figure 8.109 — DNA microarray hybridized to Panomer 9 random oligodeoxynucleotides

Figure 8.110 — Microarray spotted with the Alexa Fluor 546 and Alexa Fluor 647 Panomer 9 random oligodeoxyribonucleotides

Figure 8.111 — Secure-Seal gasket for use with microarray hybridizations

Figure 8.112 — HybriSlip covers for hybridization

Figure 8.113 — Schematic representation of molecular beacons

Figure 8.114 — Schematic representation of wavelength-shifting molecular beacons

Figure 8.115 — Schematic representation of real-time PCR with TaqMan primers

Figure 8.116 — Schematic representation of real-time PCR with Scorpion primers

Figure 8.117 — Schematic representation of real-time PCR with UniPrimers

Figure 8.118 — Four-color image of mouse intestine cryosection.

Figure 8.119 — SYTOX Blue/DNA

Figure 8.120 — Sexual reproduction of *Tetrahymena thermophila*. Texas Red(R)-X goat anti-mouse IgG (H+L) conjugate and SYTOX(R) Green nucleic acid stain.

Figure 8.121 — Macrophages treated with interleukin-4. Rhodamine phalloidin and YO-PRO(R)-1 iodide.

Figure 8.122 — Bovine pulmonary artery endothelial cells. SYTOX(R) Green nucleic acid stain and BODIPY(R) TR-X phalloidin.

Figure 8.123 — Human cheek epithelial cells. Alexa Fluor(R) 350 wheat germ agglutinin and SYTOX(R) Green nucleic acid stain.

Figure 8.124 — Nuclear deformation of an apoptotic cell visualized with the SYTOX(R) Green dye.

Figure 8.125 — Mouse brain section stained with NeuroTrace(R) 435/455 blue-fluorescent Nissl stain and counterstained with nuclear yellow.

Figure 8.126 — Binucleate bovine pulmonary artery endothelial cell. Anti- $\alpha$ -tubulin, biotin-XX conjugate, Alexa Fluor(R) 568 streptavidin and nuclear yellow.

Figure 8.127 — T1323; true blue chloride

Figure 8.128 — Zebrafish retina. TSA Kit #2, Alexa Fluor(R) 350 wheat germ agglutinin conjugate and TOTO(R)-3 nucleic acid stain.

Figure 8.129 — TOTO-3/DNA

Figure 8.130 — TO-PRO-3/DNA

Figure 8.131 — *Salmonella* gill filaments. The SYTOX(R) Green nucleic acid stain, propidium iodide or DAPI, which are all components of the FISH Counterstain Kit #2.

Figure 8.132 — Freshwater snail embryo. SYTO(R) Blue fluorescent nucleic acid stain, a component of the SYTO(R) Blue Fluorescent Stain Sampler Kit.

Figure 8.133 — Polymorphonuclear leukocytes. SYTO(R) 13 live-cell nucleic acid stain.

Figure 8.134 — Chromatin in PtK2 cells. DAPI.

Figure 8.135 — *Drosophila* embryos. Lissamine rhodamine B sulfonyl chloride and SYTOX(R) Green nucleic acid stain.

Figure 8.136 — Bovine pulmonary artery endothelial cells. Anti- $\alpha$ -tubulin, biotin-XX goat anti-mouse IgG antibody, Marina Blue(R) streptavidin and TO-PRO(R)-3.

Figure 8.137 — Mitotic spindles from sea urchin eggs. YOYO(R)-1 iodide and Texas Red(R) goat anti-rat IgG (H+L) antibody.

Figure 8.138 — Using the YOYO-1 dye to follow cell division in a sea urchin egg

Figure 8.139 — Centromeres in human metaphase chromosome spread. Fluorescein streptavidin, biotinylated anti-fluorescein and propidium iodide.

Figure 8.140 — Human metaphase chromosomes. DAPI.

Figure 8.141 — Human metaphase chromosomes. Chromosome Banding Kit #2 'SYTOX(R) Green nucleic acid stain/Methyl Green'.

Figure 8.142 — Mouse brain section stained with NeuroTrace(R) 515/535 yellow-fluorescent Nissl stain and counterstained with nuclear yellow.

Figure 8.143 — Mouse cerebellum section. NeuroTrace(R) 530/615 red-fluorescent Nissl stain and nuclear yellow.

Figure 8.144 — T2 phage genomic DNA. YOYO(R)-1 nucleic acid stain.

Figure 8.145 — Electrophoretic bandshift assays using SYBR Green I nucleic acid gel stain

Figure 8.146 — Titration of *lac* operator DNA with *lac* repressor protein.

Figure 8.147 — HL-60 cells in a comet assay. SYBR(R) Green I nucleic acid gel stain.

Figure 8.148 — Aldehyde-reactive probe (ARP) used to detect DNA damage

Figure 8.149 — Schematic diagram of the mechanism used in the EnzChek Reverse Transcriptase Assay Kit

Figure 8.150 — Detection of HIV reverse transcription using the EnzChek Reverse Transcriptase Assay Kit


## List of Technical Notes and Product Highlights

Note 8.1 — Product Highlight: Alexa Fluor Dyes for Labeling Nucleic Acids

Note 8.2 — Product Highlight: SYBR Green I and SYBR Green II Dyes for Quantitating Nucleic Acids on Plastic Wrap or Paraffin Sheets

Note 8.3 — Product Highlight: SYBR Safe DNA Gel Stain

Note 8.4 — Product Highlight: Procedure for Staining Microarrays with SYBR Green II Dye

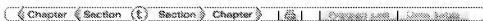

[Site Map](#) | [Change Country](#) | [Contact Us](#) | [Privacy Policy](#) | [Terms & Conditions](#) | [Online Price Policy](#) | [Feedback](#) | [Careers](#)

Copyright © 2008 Invitrogen Corporation. All rights reserved.